prorelax[®] BodyShaper

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prorelax[®]

Distributed in Europe by:
EUROMEDICS GmbH
D-53343 Wachtberg, Germany
www.euromedics.de
info@euromedics.de

prorelax[®] BodyShaper

OPERATION MANUAL

- Sustainable tightening of the skin
- Intensive massage against Cellulite
- A visible more beautiful appearance of the skin
- Visible better figure

Dear Customer,

Thank you for purchasing our **prorelax[®]direct**.

This product is designed for people to tighten their skin and bring her figure in shape. We want you to be totally satisfied with this product.

In order to achieve all the benefits from the prorelax[®]BodyShaper, we ask you to carefully read the following safety information and warnings.

Table of Contents

Table of Contents	
Intended Use & Features.....	14
Explication of the use of the attachments.....	15
Use.....	16
Warning.....	16
Technical Data.....	17
Maintenance.....	18
• Cleaning	
• Storage	
Guarantee provision.....	19
Contact.....	20

prorelax® **BodyShaper**

Intended Use


Full-body massage for a better shape and beautiful skin

- Profounding massage to tighten the skin and strengthen the connective tissue
- Easy in handling, even at problematic zones
- Especially suitable for the treatment of Cellulite

Features

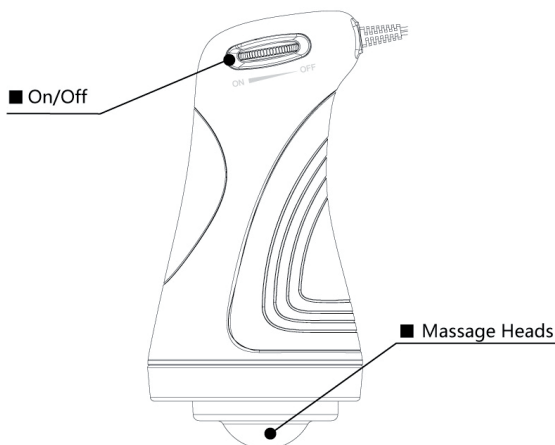- Sustainable effect because of powerful rotation of the massage head
- Extremely quiet, high tork engine, the intensity is variable
- Innovative and ergonomic form, easy in handling
- Four different formed attachments for different types of massage (see below)
- High quality, two years of full guarantee

Explication of the use of the attachments

 <p style="text-align: center;">1.</p>	 <p style="text-align: center;">2.</p>	 <p style="text-align: center;">3.</p>	 <p style="text-align: center;">4.</p>
<p>1. Ball massage head:</p>	<p>for the massage of the neck and to relax the muscles</p>		
<p>2. Double ball massage head:</p>	<p>intensive massage of the tissue</p>		
<p>3. Massage head with four balls:</p>	<p>especially suitable to improve the blood circulation</p>		
<p>4. Tai Chi massage head:</p>	<p>especially suitable for oil massages with deep action</p>		

- The protective cover is to prevent the penetration of hair into the device.
- The protective cover is not for the use with ball massage heads.

Use as shown:

Warning

- Do not use the device during pregnancy, on wounds, inflamed skin, diseased skin, at heart diseases etc.
- Wait one hour after meals.
- Children or disabled persons are to be supervised from an adult while using the device.
- Do not expose the device neither to high temperatures nor high humidity.
- Do not use the device in humidity rooms as for example in the bathroom.
- Pull the main plug after the use.
- Do not pull on the power cable of the device.
- Do not let water, oil or other liquids enter the device. (The massage head can be cleaned with water.)
- Do not use acid or metalliferous cleaning agents for cleaning the device.
- This device is not for the use of persons with reduced physic, sensory or mental ability or with a lack of experience or knowledge (including children), if there isn't a supervision or instruction regarding the use of this device by a person which is responsible for the safety of this person.
- Children are to be supervised to secure that they are not playing with the device.
- If the main plug is damaged, it has to be removed by the manufacturer, service department or a similar qualified person, to avoid danger.

Technical Data

Product name:	prorelax Bodyshaper
Model no.:	MC0313
Input power:	230 V 50-60 Hz
Power consumption:	55 W

Maintenance

Cleaning:

Be sure that the device is disconnected from the power supply and let it cool down before the cleaning. Use a smooth, lightly moistened cloth to wipe it. The device is never allowed to be in contact with water or other liquids.

Storage:

Pull the main plug out of the power socket and let it cool down, before you keep it in its box on a clean and dry place. Do not coil the cable around the device, because this may lead to wear and breakage of the cable. Do not hang it up on its cable.

GUARANTEE PROVISION

We guarantee completeness and no defects at the time of the selling.

Repairs can be done only by importers, service centres or specialist suppliers which are authorised by the EUROMEDICS GmbH.

The time of guarantee is two years from the date of the buying under the condition that the use occurs as described in the operating manual. Parts which are subjected to normal tear (for example batteries) are not included in this guarantee.

Defects within the time of guarantee have to be told to the specialist supplier or rather the service centre, which will be determining the cause. If a repairing is not possible, the product will be replaced.

The duration of the repairing is depended on the type of the defect.

The costs for delivery and collection of the device have to be paid by the buyer.

The guarantee which has to be paid by the EUROMEDICS GmbH is restricted to the amount of the purchase price. Any claims of the Buyer over and above this are excluded.

Warranty claims become void under the following circumstances:

- a) if the proof of purchase is not available,
- b) the buyer or another not authorised person has tried to repair the device,
- c) the defect was caused by the improper use by the buyer, especially by not regarding the operating manual,
- d) at obvious damage by fluids or dropping,
- e) violence while using the product.

prorelax[®]

Distributed in Europe by:
PROMEDICS HEALTHCARE GmbH
D-53343 Wachtberg, Germany
www.euromedics.de
info@euromedics.de

Article number: 705105395531